



913 pH/DO Meter

914 pH/DO/Conductometer

Combined measuring instruments for dissolved oxygen, pH, and conductivity



Handheld excellence for your most important quality parameters

The latest generation of our proven range of meters has become more powerful. Using the 913 pH/DO Meter or the 914 pH/DO/Conductometer, you can determine the pH value, concentration of dissolved oxygen, and conductivity fast, easily, and reliably. Count on trusted Swiss quality and enjoy the ease of use of our latest handheld meters.



GLP-compliant documentation

Print out your measuring results directly in the lab, or save them as a CSV report.

- Flexible: the instrument supports thermoprinters, A4 printers, or impact printers
- Simple: CSV reports for data management

Measure pH and dissolved oxygen in parallel

With the **913 pH/DO Meter**, you can measure the pH value and the concentration of dissolved oxygen simultaneously in the same sample.


- Reliable: the instrument informs you when the sensor cap needs to be exchanged
- Flexible: ideal for measurements in the field or in the laboratory



Optical measurement of dissolved oxygen

The new **O₂-Lumitrode** sensor provides faster results for dissolved oxygen than conventional sensors.

- Fast: the O₂-Lumitrode is the sensor for dissolved oxygen with the fastest response time on the market
- Maintenance-free: the O₂-Lumitrode sensor is always ready to use
- Robust: the sensor membrane is protected from mechanical stress and interfering external light

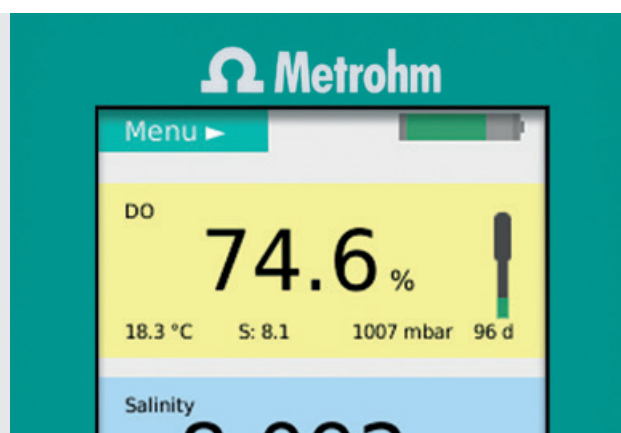

Measure pH value, dissolved oxygen, and conductivity with the same instrument

With the **914 pH/DO/Conductometer**, you can measure all three parameters (pH, dissolved oxygen, and conductivity) on the same platform.

- Accurate: when measuring dissolved oxygen in solutions containing high salt concentrations, this can be compensated automatically
- Dependable: the instrument informs you when the sensor cap needs to be exchanged

Reliable measuring results

The instrument informs you when sensor performance declines and the sensor cap needs to be replaced.

- Convenient: active monitoring of sensor performance
- Cost-optimized: there is no defined limit to the lifetime of the sensor. Instead, it depends on the intensity of usage
- Reliable: correct results over the entire lifetime of the sensor cap

ORDERING INFORMATION

913 pH/DO Meter	
29130020	pH/DO Meter
29130120	pH/DO Meter with Primatrode and O ₂ -Lumitrode in a case with various accessories
29130220	pH/DO Meter with stand
914 pH/DO/Conductometer	
29140030	pH/DO/Conductometer
29140130	pH/DO/Conductometer with conductivity sensor and O ₂ -Lumitrode in a carrying case with various accessories
29140230	pH/DO/Conductometer with stand
O ₂ -Lumitrode	
61116000	O ₂ -Lumitrode with calibration vessel and oxygen standard 0 %
62326000	Oxygen standard 0 %
65623000	Replacement set for O ₂ -Lumitrode with calibration vessel and oxygen standard 0 %